

**A LARGE
EMERGING
MARKET**

**53
MILLION
Target
Population
in US alone**



INVEST in the rapidly
growing VAGINAL
HEALTH market

EXPANDABILITY

diVa® operates on the JOULE™ platform, the highest quality and most expandable energy based system presently on the market. When you own the JOULE system, you have invested in the future of your practice. The JOULE offers the opportunity to, overtime, add new revenue streams with over a dozen in-demand applications in one system. It not only helps to grow your practice, but JOULE grows with your practice.

INNOVATIVE AND RELIABLE

Sciton, Inc. is an employee owned company, focused on quality and innovation. 95% of all Sciton systems ever sold since 1997 are still in service.

ONLY SCITON HAS A COMPLETE SOLUTION

A complete partnership, a novel concept committed to your success.

POST SALE SUPPORT

Sciton prides itself on partnering with you from the beginning and beyond.

PRACTICE DEVELOPMENT PARTNER (PDP)

A live dedicated person to help your practice.

- Ongoing research focus in best practices, marketing pearls, success stories, and competitive structure
- Design and implement creative and innovative ideas supporting practices to help you stay on top
- Partner closely with you to provide solutions to obstacle
- Act as liaison to entire staff, MD, RN, PA, & Front Desk



SCITONPRO.COM

Sciton professional website provided free of charge:

- Digital practice marketing assets
- Available 24/7

SMART START™ GUIDE

5 steps to operating a thriving aesthetic practice:


- STEP 1: Get launch ready
- STEP 2: Become an expert
- STEP 3: Prepare your staff and office
- STEP 4: Know your patient
- STEP 5: Develop your marketing plan



SUCCESS BUILDER™

Access to experienced third-party professionals and tools to support your practice growth.

- Clinical education
- Marketing tools
- Business development
- Web or video assets

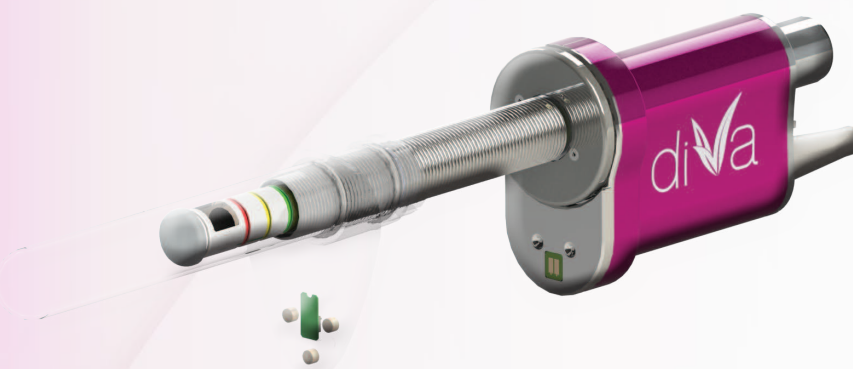

diVa®

diVa® an innovation in
laser resurfacing technology
for the improvement of
gynecological and
genitourinary/urology health



925 Commercial Street, Palo Alto, California 94303
Phone: (888) 646-6999 • Email: info@sciton.com

www.sciton.com • f in

©2019 Sciton, Inc. All rights reserved. Sciton and diVa are registered trademarks, JOULE and HALO are trademarks of Sciton, Inc.
Printed in USA 2600-029-23 Rev. D

The JOULE 2940 nm fractional delivery system has been cleared for skin resurfacing, and ablation and coagulation of soft tissue. Genitourinary surgery of soft tissue, including: treatment of bladder, urethral and ureteral tumors; superficial urinary bladder tumors; invasive bladder carcinomas; urethral and penile hemangioma; urethral strictures; lesions of the external genitalia; condylomas; bladder neck obstructions. Gynecology indications include cervical intraepithelial neoplasia (CIN), herpes simplex, endometrial adhesions, cysts and condylomas. The JOULE1470 nm fractional delivery system has been FDA cleared for skin resurfacing and coagulation of soft tissue.





diVa®, the LATEST INNOVATION developed by Sciton, EXPANDS the BROAD RANGE OF THE JOULE™ platform to address woman's vaginal health.

DIVA TECHNOLOGY

HYBRID FRACTIONAL LASER (HFL) TECHNOLOGY

diVa® utilizes the revolutionary Hybrid Fractional Laser technology developed for HALO™, the world's first hybrid fractional laser, to provide tunable independent levels of ablation and coagulation, as well as varying levels of density. This allows the provider to customize settings based on the patient's symptoms.

HIGH PRECISION AUTOMATION (HPA)

The motorized guidance system allows providers to deliver energy with increased consistency, limited user variability, and decreased procedure time. HPA technology promotes treatments that are typically uniform, predictable, and fast 3-5 minutes.

PATIENT COMFORT

The single-use Strengthened Quartz Dilator (SQD) consists of a quartz tip that remains stationary to expand and flatten the vaginal wall tissue for a uniform and typically comfortable treatment – no manipulation of sensitive tissue, the HPA rotates within the SQD. As the handpiece delivers energy to the treatment area, pulses are delivered through a new section of quartz to promote laser efficiency. The disposable SQD minimizes the risk of cross contamination and shortens downtime between treatments.

CUSTOMIZE ENERGY TO THE TREATMENT AREA

The diVa hand piece allows for customized energy delivery at different angles to address specific areas of concern which allows the provider to further optimize the treatment. This gives diVa the unique ability to provide individualized energy based on the type of treatment and thickness of the vaginal wall.

WAVELENGTHS THAT ARE PURPOSE DRIVEN

diVa delivers both ablative and non-ablative wavelengths to the same microscopic treatment zone: 2940 nm for ablation and 1470 nm for coagulation. The 2940 nm wavelength is used to deliver 100 to 800 microns depth of ablation and the 1470 nm wavelength is used for delivering 200 to 700 microns depth of coagulation to the vaginal soft tissue. Choose each wavelength's parameters independently to create a unique synergistic effect.



EASY TOUCH-SCREEN TREATMENT USER INTERFACE

DIVA IS INDICATED FOR THE TREATMENT OF SCARS

If your patients have experienced vaginal childbirth or any gynecological surgery, chances are, they have internal and/or external vaginal scarring and issues associated with those scars. Sciton's suite of women's health options, including diVa, will end the discomfort associated with those scars.




Courtesy of Michael Coyle, DO



Courtesy of Michael Coyle, DO

STRENGTHENED QUARTZ DILATOR (SQD)

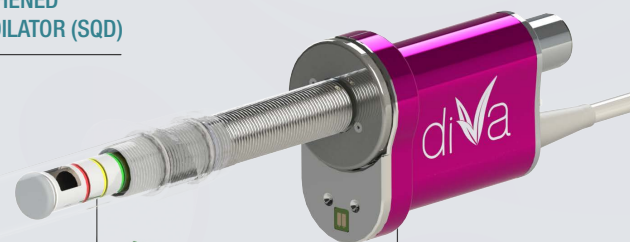

DELIVERY AT MULTIPLE ANGLES

HIGH PRECISION AUTOMATION (HPA)